

ELISA Flex: Hamster IgG (HRP)

3154-1H-6 |

ELISA Flex kit for quantitative determination of native Syrian hamster IgG (*Mesocricetus auratus*, golden hamster) in solution, e.g. serum/plasma samples or cell supernatants.

| The kit includes | | 3154-1H-6 for 6 plates | |
|--|-----------------------------|-----------------------------------|--|
| Capture mAb: | MTAE142 (0.5 mg/ml) | 300 µl | |
| Detection mAb: | MTAE913, biotin (0.5 mg/ml) | 150 µl | |
| Streptavidin-HRP | | 80 µl | |
| Native Syrian hamster IgG ELISA standard | | 1 vial | |
| Standard reconstitution buffer A5 | | 1 ml | |

To ensure total recovery of the stated quantity, vials have been overfilled.

Shipping and storage

Shipped at ambient temperature. All reagents should be stored at 4-8 °C upon receipt, except the standard which should be stored at -20 °C. Antibodies are supplied in sterile-filtered PBS with sodium azide (0.02%). Streptavidin-HRP is supplied in PBS with 0.002% Kathon CG. The expiry date indicates how long unopened products, stored according to instructions, are recommended for use.

General and Preparations

Specificity

The kit contains a matched pair of monoclonal antibodies (mAbs) specific for Syrian hamster IgG (*Mesocricetus auratus* golden hamster). No cross-reactivity with IgG from Armenian hamster.

Standard range

1-100 ng/ml

Calibration


No international standard exists for calibration.

Reconstitution of ELISA standard

Reconstitute the ELISA standard to a stock solution of 10 µg/ml by adding 1 ml of the standard reconstitution buffer. Allow the standard to dissolve for 5 minutes and mix thoroughly. The standard should be kept in aliquots at -20 °C. Avoid repeated freeze-thaw cycles.

Preparation of standard curve

Prepare within 30 minutes of use. Volumes are sufficient for duplicates.

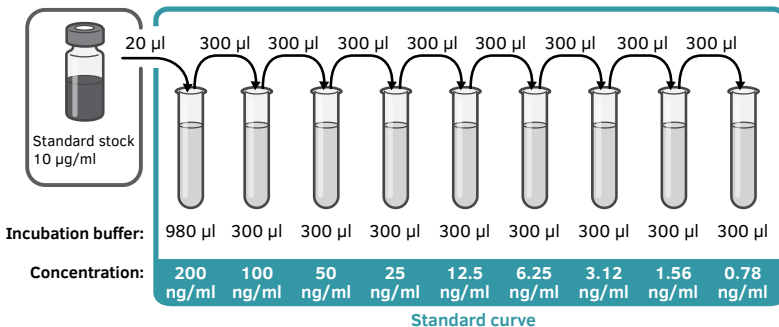


Protocol

Day 1

1. Add 100 μl /well of capture mAb MTAE142 diluted to 2 $\mu\text{g}/\text{ml}$ in PBS, pH 7.4. Use high protein binding ELISA plates. Incubate overnight at 4-8 $^{\circ}\text{C}$.

Day 2

2. Empty the plate and add 200 μl /well of PBS with 0.05% Tween 20 and 0.1% BSA (incubation buffer) to block the plate. Incubate for 1 hour at room temperature.
3. Wash the plate 5 times with PBS containing 0.05% Tween 20 (300 μl /well).
4. Add 100 μl /well of samples or standards diluted in incubation buffer. Include assay background control, i.e. wells without standard. Incubate for 2 hours at room temperature.
5. Wash as above.
6. Add 100 μl /well of detection mAb MTAE913-biotin diluted to 1 $\mu\text{g}/\text{ml}$ in incubation buffer. Incubate for 1 hour at room temperature.
7. Wash as above.
8. Add 100 μl /well of Streptavidin-HRP diluted 1:1000 in incubation buffer. Incubate for 1 hour at room temperature. Please note that sodium azide used in buffers will inhibit HRP activity.
9. Wash as above.
10. Add 100 μl /well of TMB substrate (product code: 3652-F10) and incubate at room temperature, protected from direct light for 15 minutes.
11. Add 100 μl /well of 0.2 M H_2SO_4 to stop the reaction.
12. Measure the optical density in an ELISA reader at 450 nm within 15 min. Preferably use a reader capable of subtracting a reference wavelength of between 570 and 650 nm. Representative standard curve shown below.

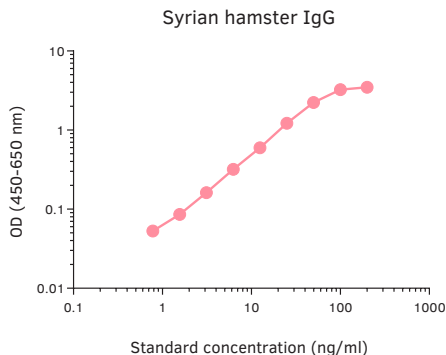

Developed and manufactured by MABTECH AB, Sweden, whose quality management system complies with the standards ISO 9001:2015 & ISO 13485:2016.



The products are for research use only.

MABTECH shall not be liable for the use or handling of the product or for consequential, special, indirect or incidental damages there from.

Mabtech AB (Head Office)

Sweden
Tel: +46 8 716 27 00

mabtech@mabtech.com

Mabtech, Inc.
USA

www.mabtech.com