

# ELISA Flex: Human IL-27 (ALP)

3458-1A-6 | 3458-1A-20

ELISA Flex kit for quantitative determination of native and recombinant human IL-27 in solution, e.g. cell supernatant.

<b>The kit includes</b>	<b>3458-1A-6 for 6 plates</b>	<b>3458-1A-20 for 20 plates</b>
Capture mAbs: MT27/361 (0.5 mg/ml)	300 µl	1000 µl
Detection mAb: MT140, biotin (0.5 mg/ml)	80 µl	250 µl
Streptavidin-ALP	80 µl	250 µl
Recombinant human IL-27 ELISA standard	1 vial	1 vial
Standard reconstitution buffer A5	1 ml	1 ml

To ensure total recovery of the stated quantity, vials have been overfilled.

## Shipping and storage

Shipped at ambient temperature. All reagents should be stored at 4-8 °C upon receipt, except the standard which should be stored at -20 °C. Antibodies are supplied in sterile-filtered PBS with sodium azide (0.02%). Streptavidin-ALP is supplied in 0.1 M Tris buffer with 0.002% Kathon CG. The expiry date indicates how long unopened products, stored according to instructions, are recommended for use.

# General and Preparations

## Specificity

The kit contains a matched pair of monoclonal antibodies (mAbs) specific for native and recombinant human IL-27.

## Standard range

32-3200 pg/ml

## Calibration

No international standard exists for calibration.

## Analysis of serum and plasma samples


Analysis of serum/plasma requires the use of ELISA diluent (product code: 3652-D2). The ELISA diluent blocks heterophilic antibodies, commonly found in serum/plasma, from cross-linking the assay antibodies, thereby preventing false positive read-outs. The ELISA diluent should be used for dilution of standard, samples, and detection antibody.

## Reconstitution of ELISA standard

Reconstitute the ELISA standard to a stock solution of 0.5 µg/ml by adding 1 ml of the standard reconstitution buffer. Allow the standard to dissolve for 5 minutes and mix thoroughly. The standard should be kept in aliquots at -20 °C. Avoid repeated freeze-thaw cycles.

## Preparation of standard curve

Prepare within 30 minutes of use. Volumes are sufficient for duplicates.

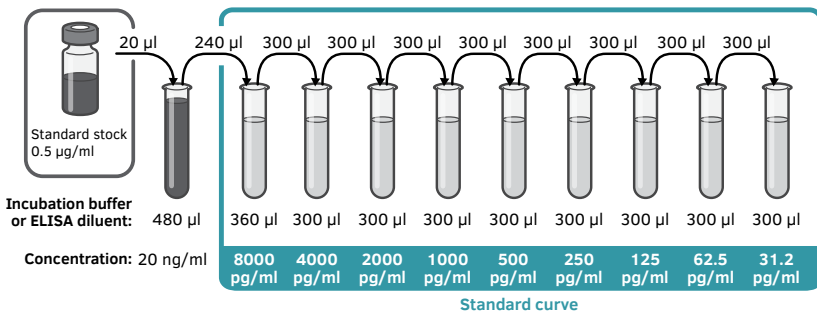


# Protocol

## Day 1

1. Add 100  $\mu$ l/well of capture mAbs MT27/361 diluted to 2  $\mu$ g/ml in PBS, pH 7.4. Use high protein binding ELISA plates. Incubate overnight at 4-8  $^{\circ}$ C.

## Day 2

2. Empty the plate and add 200  $\mu$ l/well of PBS with 0.05% Tween 20 and 0.1% BSA (incubation buffer) to block the plate. Incubate for 1 hour at room temperature.
3. Wash the plate 5 times with PBS containing 0.05% Tween 20 (300  $\mu$ l/well).
4. Add 100  $\mu$ l/well of samples or standards diluted in incubation buffer or ELISA diluent. Include assay background control, i.e. wells without standard. Incubate for 2 hours at room temperature.
5. Wash as above.
6. Add 100  $\mu$ l/well of detection mAb MT140-biotin diluted to 0.5  $\mu$ g/ml in incubation buffer or ELISA diluent. Incubate for 1 hour at room temperature.
7. Wash as above.
8. Add 100  $\mu$ l/well of Streptavidin-ALP diluted 1:1000 in incubation buffer. Incubate for 1 hour at room temperature.
9. Wash as above.
10. Add 100  $\mu$ l/well of pNPP substrate (product code: 3652-P10) and incubate at room temperature protected from direct light for approximately 60 minutes.
11. Measure the optical density in an ELISA reader at 405 nm. Preferably use a reader capable of subtracting a reference wavelength of between 570 and 650 nm. Representative standard curve shown below.

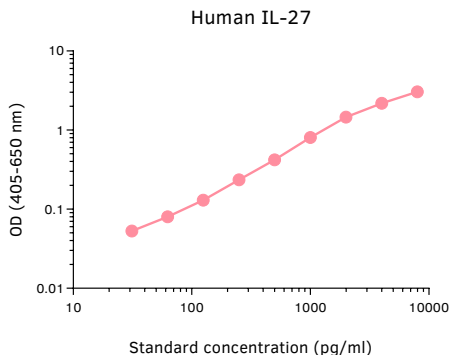

Developed and manufactured by MABTECH AB, Sweden, whose quality management system complies with the standards ISO 9001:2015 & ISO 13485:2016.



**The products are for research use only.**

MABTECH shall not be liable for the use or handling of the product or for consequential, special, indirect or incidental damages there from.

---

**Mabtech AB (Head Office)**  
Sweden  
Tel: +46 8 716 27 00  
mabtech@mabtech.com

**Mabtech, Inc.**  
USA  
Tel: +1 513 871-4500  
mabtech.usa@mabtech.com